



 Products & Capabilities

Alina[®]

User-friendly disposable pen injector platform
for variable and multi-dose treatments

 Stevanato Group

MANAGING COMPLEXITY, DELIVERING VALUE



Alina 60
Pen injector Platform

3mL

0

About Alina[®]

Alina[®] is an easy-to-use disposable pen injector platform with variable dosing.

The pen injector's performance, features and intuitive handling are compatible with established therapeutic regimens as well as innovator drug therapies that require a multidose delivery device such as diabetes and weight management.



Features

Suitable for **ISO 3mL glass cartridges** & optimized for Stevanato Group's range of glass solutions

Compatible with **standard** pen needles

Dose indicator, with registered design, for **improved readability**

Removable cap

Dial-back **dose error correction**



Benefits



Ideally suited for
diabetes and weight
management



Visual and audible
feedback for dose setting,
correction & injection



Tested & verified
according to ISO11608
standards




Customization available
for primary packaging &
therapy needs

Handling Steps




1 Remove cap

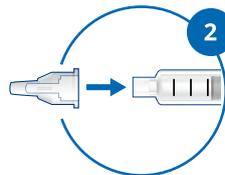

2 Attach the needle



3 Hold the device
upright and prime

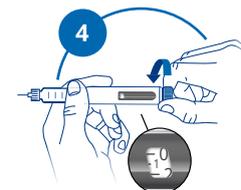

4 Turn the dial
to select the dose



5 Inject carefully
into a suitable
injection pad




6 Remove the needle
& recap the device

Your Integrated Solutions Provider

Production of Alina® takes place at our state-of-the-art, ISO13485 certified, FDA-inspected manufacturing site in Germany.

We have invested in component tooling, automated sub-assembly equipment and optimized production processes to accelerate product launches and minimize up-front investments for customers.

Additionally, we have a range of optional products and services to further simplify your supply chain.

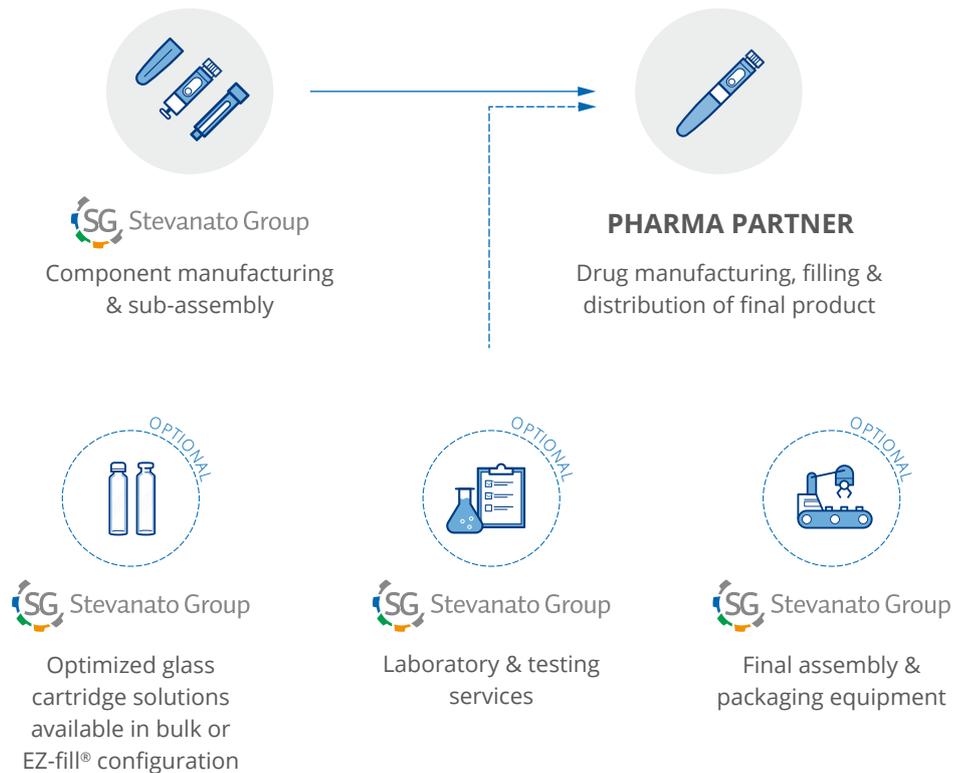



Headquarters

Via Molinella, 17
35017 Piombino Dese
Padova, Italy

**Global presence with operating units
and commercial offices in 10 countries.**

Visit our website to contact us.



SG, Stevanato, Vertiva, EZ-fill, Nexa, Alba, Alina, Fina are registered trademarks of Stevanato Group and/or its affiliates in one or more jurisdictions.

stevanatogroup.com

ALINA_202412